

Multi-function scraping detoxifying body massager

Scraping artifact



Due to changes in product update models,
the color function is based on actual products.

PRODUCT PRINCIPLE

This product is a collection of negative pressure suction, magnetic wave therapy, fever and blood pressure, acupressure, EMS myoelectric stimulation, far infrared therapy and other functions. It can be used for scraping and sucking, discharging acid and dehumidifying, relaxing nerves, passing through vitality, adjusting function, dispelling cold and dehumidifying, eliminating fatigue, clearing fat and losing weight, tightening skin and so on. It is a new generation of health equipment for every family.

Pulsating energy therapy can quickly clear the lymphatic vessels and meridians of the body, and thus play the role of health and weight loss beauty. The acupoint massage of the bump forms a small squeeze and eliminates fat, which dissipates hard fat and is very effective in eliminating cellulite. The use of one-stop communication, regulation, and health care, all-round treatment to improve various sub-health.

The instrument has an aerobic effect on the fat in certain areas of the body, stimulating the body's meridians and acupoints, and simultaneously activating the cell viability in the region, which can improve the metabolism of the whole body; the heat energy can raise the temperature of the deep fat to make it Dissolves metabolism and excretes excess toxins and moisture, tightening the collagen fibers of the dermal layer of the skin. While lowering the level of fat, it can also tighten the skin and shape the body curve.

NEGATIVE PRESSURE GRAVITATION:

- A. Aerobic exercise for upward absorption of skin and adipose tissue
- B. Enhance lymphatic detoxification' function and blood circulation
- C. Elimination of toxins through the lymphatic circulation.
- D. Scraping.
- E. Cupping

PHYSICAL HYPERTHERMIA:

- A. Proliferate dermal collagen by thermal effects.
- B. Continuous heat, nerve relaxation, pressure release

MAGNETIC WAVE THERAPY:


- A. Promote inflammation and eliminate inflammation, swelling and pain
- B. Eliminate fatigue and promote physical recovery.

EMS MICRO-ELECTRIC STIMULATION:

- A: Stimulate the acupuncture points by micro-current, so that the meridians can be opened in an instant to achieve the purpose of passing through the active.
- B: The micro-current stimulates the peripheral sensory nerve, shrinks the pores, and the skin is delicate and lustrous.

FAR INFRARED THERAPY:

- A. Enhance the body's immunity and the ability of tissue regeneration of biological cells to accelerate the supply of nutrients and enzymes to promote good health.
- B. It can stimulate and promote the metabolism of cells. Elimination of venous congestion improves blood circulation, enhances cell activity and stimulates collagen regeneration to make your skin smooth and tender.

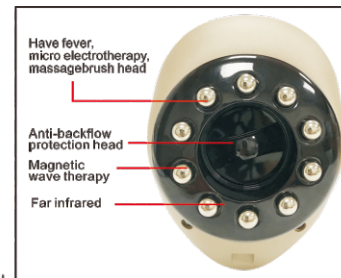

Have fever, micro electrotherapy, massage brush head

Anti-backflow protection head

Magnetic wave therapy

Far infrared

HOW TO USE

Technical specifications
Input voltage: 100V-240V
Power: 20W
Working voltage: 9V
Packing size: 262x152 x 110
Material body: abs, 304 stainless steel
Made in china

Control panel introduction



Biological micro-electric button

Power on/off button

Negative pressure start button

Hot button

SCRAPER OPERATION INSTRUCTIONS

Functional module	Icon display	Scraper
Switch button	/	Click on the boot and then press the shutdown
EMS micro-electric	EMS	1. EMS has 0-4 files, the corresponding LCD screen displays 1-4 grids, and adjusts the gear position according to EMS cycle 0-4 files-4-0 files.
Negative pressure suction	VACUUM	1. VAC has 0-4 files, the corresponding LCD screen displays 1-4 grids, press the negative pressure button cycle 0-4 files-4-0 gears to adjust the gear position
Heating	Heat	1. Heating switch: press to turn on the heating, heat the icon display, then press the lower heating is off, the heating icon is off 2. Heating setting: press and hold the heating button for 1.2 seconds to enter the temperature setting, the LCD screen temperature display flashes, the adjustment range is 35-50 degrees Celsius, press the heating button to adjust the temperature cyclically, adjust the temperature by 1 degree Celsius, and adjust the temperature. the temperature shows the corresponding temperature setting and flashes, 25 does not operate save to exit the setting state, the temperature display returns to the correct temperature 3. During heating starting or during the process, do not heat when the temperature is not displayed.
Backlight	/	TOSS does not operate, the backlight automatically turns off, touches the button backlight to turn on (function does not operate)
Charging	Battery	1. Display the lightning bolt icon when charging, the power stall will display the current power in real time, and work when charging.
Red light	Light	The red light is on when the machine is working, and the red light is off when not working.
Automatic shut-down	/	No operation automatically shuts down within 30 minutes after power on.

COMMON PROBLEM

When the following occurs during use, please refer to the corresponding processing method below for processing.

1. **The power supply is unreachable**
 1. Do not press and hold the power button, please press and hold this button for more than 1 second.
 2. If it is a charging model, please charge it in about 3 hours if it is not enough.
 3. power automatically turns off the machine does not use automatic shutdown 30 minutes after turning on
2. **Can't boot work**
 - Poor contact between the power adapter and the instrument socket
3. **Micro current is not available**
 1. If the part is not in full contact with the skin or the skin is dry, please contact the skin in the main operation head to keep the skin moist. Apply cosmetic cream and gel-like products to keep the skin moist.
 2. micro-current energy is different for each person's skin feeling, please adjust the energy according to personal strength
 3. Sometimes feel the skin tingling, please adjust the energy level to keep the skin moist and fully contact the skin.

CHARGING INSTRUCTIONS (CHARGE ONLY)

If it is the first time you use the charging, you must charge it for 3 hours. The charging time is about 70 minutes. Depending on the mode selected, the usage environment will change in time. This product is equipped with security protection and will automatically shut down after 30 minutes of power on. When charging, there will be a lightning bolt icon in the upper right corner of the screen, and the power will be automatically turned off after the machine is fully charged.

Finally explain the main:
Different versions have incomplete functions, and the specific functions are subject to the actual style.

ABOUT MAINTENANCE

After each use, the surface of the machine must be cleaned. When cleaning, you must turn off the power and unplug the power plug. If you encounter stubborn stains, you can wipe it with a sponge or cloth moistened with a diluted neutral detergent. Please do not use strong alkalis. Volatile oil, delustering liquid, strong essential oil, medicated oil, etc. can be wiped locally with a small amount of alcohol. It is recommended to use with the warm wormwood essential oil provided by our factory. After cleaning, please put on the bottom cover. This instrument is not waterproof. You cannot use it in water or wash the host and power supply with water.

Note: Due to the irregular upgrade of the product, the model number of the manual may be different from the actual product. This manual is for reference only.

After cleaning, remove the water-absorbing hose, turn on the suction and run for 5-10 minutes, and drain the water without noise. If the cleaned product makes abnormal noise during use, it is because the water in the pipe and pump has not been drained, just turn on the suction and drain it without load. Even if it is not drained, it does not affect performance and normal use.

INSTRUMENT INSTRUCTIONS

1. Negative Pressure

For scraping and cupping, scraping: 1-2 gears, cupping: 3-4 gears. Suggestions for operation:
1. Scrape: apply essential oil, direction: from top to bottom, strength: apply evenly, heat up/out of Sha. 2. Cupping: Duration of each stay: generally about half a minute, the length of time can be adjusted according to the degree of Sha-out of the individual.

2. Microelectronics



It can replenish the body's bioelectricity, activate cells, relieve fatigue, and restore skin elasticity, so as to achieve the function of solving some skin problems and delaying skin aging. It is recommended that this function be turned on all the time, and the gear setting is determined by personal preference.

3. Heating

Setting steps: 1. Long press the heating button for 3 seconds until the temperature display flashes.
2. Click the heating button to set the temperature (35-50°C). 3. After the setting is completed, the temperature no longer flashes, that is, the setting is successful. (It is recommended to set around 40°C)

MASSAGE TECHNIQUE

The following is upward, the outward direction is the main direction, and the intensity is mainly gentle.

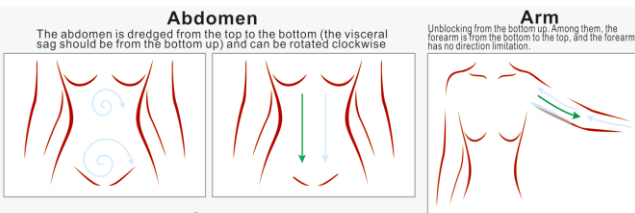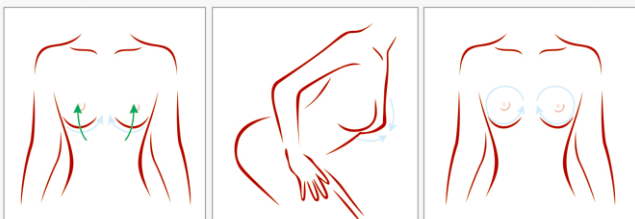

PRECAUTIONS

1. The body has varicose veins and should not be used.
2. During pregnancy or menstruation, it should be suspended.
3. Because the blood is concentrated in the digestive system after meals, it should not be used.
4. After use, the surface of the skin is dilated. It should be rested for 30-45 minutes before bathing.
5. When using, you should follow the same direction, not one on the other, one left and one right.
6. Due to the phenomenon of sputum caused by negative pressure gravitation, it should be suspended for a few days, and then use the instrument, the effect is better.


Special attention:

1. Before using the instrument, take a few drops of essential oil and apply it evenly around the skin. Massage until the skin is moisturized to avoid discomfort caused by dryness.
2. Do not apply too much essential oil to avoid product blockage.



ONE-KEY CLEANING AND EASY STERILIZATION

Upgrade the washing machine core to solve the problem of insufficient suction

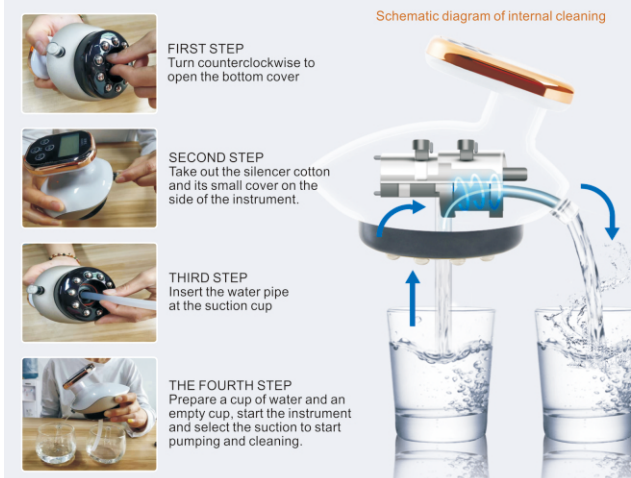

Schematic diagram of internal cleaning



FIRST STEP
Turn counterclockwise to open the bottom cover



SECOND STEP
Take out the silencer cotton and its small cover on the side of the instrument.



THIRD STEP
Insert the water pipe at the suction cup



THE FOURTH STEP
Prepare a cup of water and an empty cup, start the instrument and select the suction to start pumping and cleaning.


AFTER SALES SERVICE

Normally, when a failure occurs during use, the warranty (without consumables) is within one year from the date of purchase. If the warranty period fails, the repair is performed at the actual cost.

1. Under the premise of complying with the instructions for use of this product, use the repair if there is a failure.
 2. During the warranty period, if the following conditions are met, it is a paid repair
 - Incorrect use, incorrect repair, alteration or damage due to malfunction.
 - After purchasing the product, it falls to the fault and breaks.
 - Other natural disasters such as fire, earthquake, flood, etc. cause voltage abnormality failure and damage
 - Abnormality in connection with other electrical appliances, resulting in product failure and damage
- When items such as consumables are replaced

Warranty card

Customer name							
Distribution unit							
Contact number							
Date of purchase	Year	Month	Day	---	Year	Month	Day
Warranty period	Year	Month	Day	---	Year	Month	Day

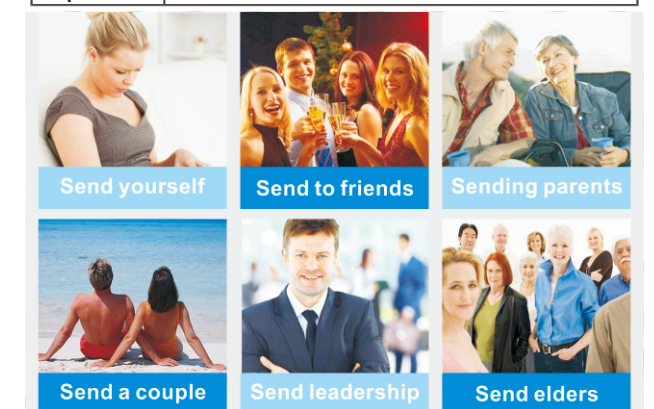

Send yourself

Send to friends

Sending parents

Send a couple

Send leadership

Send elders